

Product Catalog

2022

Table of Contents

A	D (continued)	P
All-Bite Sensor Holder..... 12	Duplicator 21	Panoramic Film 21
Aprons, Lead, lead free 14,15	Duplicating Film 21	Perfect Fit Camera Sleeves 8
Apron Hangers..... 15	E	PSP Accessories 5
B	Econo Pak X-ray Films 16-20	PSP Envelopes 2-4
Bib Apron 14,15	Econo Safe 'N' Sure 4	Perfect Pocket Mounts 24, 25
Bite Block Baggies 7, 23	Endodontic Film Holder 7, 23	R
Bite Wing Loops 5, 22	Exposure Articulating Film..... 13	RAPiD Positioner/Parts 7, 23
Bullseye Aiming Ring 12	F	S
Butterfly Barrier X-ray Film 20	Face Shields 14	Safelights 21
C	Fresh Pak X-Ray Film 16-20	Safe 'N' Sure PSP Envelopes 2-4
Cable Saver 11	G	Safe-T-First Face Shields 14
Cassettes, flexible 21	Greylight Mounts 24, 25	Sensibles Sensor Positioner..... 11
Chemistry 20	Grip-Right Mounts 24, 25	SMART PSP Positioner 6
Child's Size Aprons 14,15	H	Snuggies Sensor Sleeves..... 10
Clean Ups PSP Cleaners..... 6	Hangers for Aprons 15	Stickems Adhesive Tabs 5, 22
Comfees Sensor Sleeves 10	I	Sticky Bites 5, 22
Comfee Bytes CBCT Covers 9	Illuminators 21	SUPAs 5, 22
Comfee Ears CBCT Covers 9	K	T
Composite Dispenser Sleeve 8	Klear-Vue X-Ray Mounts 24, 25	Tranzform hanger..... 15
Cushies 5, 22	L	U
D	Lead Free Protective Aprons 26	Uni-Vers-All Positioner..... 12
Darkroom Safelights 21	Loops, Bite-Wing "Loop-EZ" 5, 22	Uni-Hanger, Aprons 15
Deluxe Cushies 5, 22	M	X
Deluxe Safe 'N' Sure 3	Mounts, X-Ray 24, 25	Xpress F Speed X-Ray Film 19-20
Digi Loops 10		X-Ray Films 16-20
Dosimeter Service..... 13		X-Ray Mounts 24, 25
Dental Illuminators 29		X-Ray Protective Apparel 14, 15

Item	Page	Item	Page	Item	Page	Item	Page	Item	Page	Item	Page	Item	Page
11532	21	18235	17,18	2316	24,25	40900	7,23	80125	8	408610	7,23	81-3456	24,25
15700	5,22	18236	17,18	23219	24,25	40901	7,23	80130	4	408620	7,23	81-4021	24,25
15710	5,22	18310	19,20	23226	24,25	40909	11	80131	3	408630	7,23	81-4040	24,25
15720	5,22	18410	19,20	23228	24,25	40910	11	80141	2	408700	6	81-4567	24,25
15725	5,22	18600	17,18	23239	24,25	40915	11	80151	2	408800	12	86011	24,25
15730	11	18601	17,18	23608	24,25	40916	12	80161	2	409150	11	81-6021	24,25
15731	11	18620	19,20	24220	20	40917	12	80171	2	409155	11	81-6040	24,25
15732	11	18621	19,20	24300	13	70104	7,23	80200	10	409260	7,23	81-6045	24,25
15733	11	18635	19,20	25200	24,25	70107	6	80201	10	409270	7,23	81-6061	24,25
15780	5,22	18636	19,20	25203	24,25	75050	15	80210	10	409280	7,23	81-6081	24,25
15800	5,22	19103	21	25608	24,25	75051	15	80211	10	409290	7,23	81-6162	24,25
15810	22	19104	21	25609	24,25	75053	15	80212	10	18095	21	81-6180	24,25
16087	5,22	19198	21	25209	24,25	75056	15	80213	10	18096	21	81-6200	24,25
16102	13	19203	21	25206	24,25	75061	15	80222	10	18095-G	21	81-7021	24,25
16220	15	19235	20	25215	24,25	75062	15	80223	10	18096-G	21	81-8142	24,25
16401	15	19236	20	25219	24,25	75090	15	80500	11	408570-S	7,23	81-0500	24,25
18060	20	19298	21	25261	24,25	80100	4	81449	9	408580-S	7,23	602-203G	24,25
18070	21	19600	20	25240	24,25	80101	3	81500	9	40909-U	11	604-209G	24,25
18071	21	19601	20	25225	24,25	80110	4	403030	7,23	40910-U	11	6045-239G	24,25
18095	21	19610	20	25226	24,25	80111	3	408500	6	409270-S	7,23	606-215G	24,25
18096	21	19620	20	25227	24,25	80120	4	408520	6	75050NL	14	608-219G	24,25
18200	17,18	19635	20	25228	24,25	80120DBX	4	408540	6	75051NL	14	610-262G	24,25
18201	17,18	19636	20	40310	7,23	80121	3	408560	7,23	75053NL	14	614-225G	24,25
18204	17,18	21502	21	40311	7,23	80121DBX	3	408570	7,23	75056NL	14	614B-625G	24,25
18206	17,18	21522	21	40312	7,23	80122	4	408580	7,23	75061NL	14	618-228G	24,25
18207	17,18	23206	24,25	40313	7,23	80123	4	408590	7,23	75062NL	14	618-230G	24,25
18215	19,20	23209	24,25	408540	7,23	80124	4	408600	7,23	75090NL	14	610-270G	24,25



Now, more than ever, this catalog is your one-stop resource for "everything imaging".

In addition to our best-selling accessories for phosphor plates, sensors, and film, we now offer a selection of products for the CBCT user.

These include our Comfee Bytes, Comfee Ear sleeves, Bite Block Baggies and CBCT lead and lead free aprons.

Also new to this catalog is our Deluxe Safe'n'Sure and Econo Safe'n'Sure in a handy wall mounted Dispenser! Speed up your psp processing by having your envelopes in a handy dispenser either on or next to your processor!

And for sensor users we have re-designed our top selling Sensibles sensor holder. It's better than ever!

All this and more is why Flow Dental has been and will be your trusted resource for everything imaging!

And if you have an idea for an imaging accessory or an idea to improve anything we already make, send me and e-mail. We would love to hear from you.

Stay safe,

William Winters

William Winters **President**

Save Money with Flow Dough!

Select Flow products contain valuable **Flow Dough Dollars!** Redeem them for free film or any one of hundreds of imaging accessories in this catalog.

Sign up and redeem your Flow Dough today at: www.FlowDough.com



Phosphor Plate Products

- PSP Envelopes
- PSP Positioning Accessories
- PSP Positioners



CBCT Products

- CBCT Bite Blocks
- CBCT Ear Sleeves
- CBCT Bite Block Baggies
- CBCT Patient Aprons



Sensor Products

- Sensor Protective Sleeves
- Sensor Positioning
- Sensor Wire Protector
- Positioner for Hand-Held X-rays



Protection Products

- Dosimeter Badge Service
- Face Shields
- Lead-free Aprons
- Lead Aprons
- Apron Hangers



X-Ray Film Products

- D & F Speed Intra-oral Films
- Barrier Packet X-ray Films
- Film Positioners
- Chemistry
- Panoramic/Duplicating Films
- X-ray Mounts
- Illuminators & Safelights



Non X-Ray Barrier Products

- Camera Sleeves
- Articulating Film
- Composite Dispenser Sleeves



Safe 'n' Sure OPT™ PSP Envelopes

Only OPT features pre-loaded cardboard inserts. By pre-loading the cardboard inserts we make loading your psp plates quick and easy. This also ensures the correct side of the plate is always facing the x-ray tube.

OPT features a semi-transparent side so you can easily verify proper plate insertion.

If you use a psp plate with magnets, the circular cut out on the cardboard insert guarantees correct plate placement. The magnet will nest neatly inside, (see the "Load" image below)

If you don't use psp plates with magnets, the cardboard insert still provides extra protection and rigidity.

And our patented EZ-Glide tab makes loading and sealing your envelopes fast and effortless. To remove your plates, our Easy-tear center seam makes for quick and easy plate removal, and without having to touch the plate during processing.

OPT comes 250 per box, and each box has your envelopes neatly stacked for easy access.



LOAD



SEAL



UNLOAD

Each envelope comes with cardboard insert already pre-loaded!

Safe 'n' Sure OPT

Size 0: 80141 Size 1: 80151 Size 2: 80161 Size 3: 80171

250/box



Whatever your preferences are, there's a Safe 'n' Sure PSP envelope that's perfect for you!

Here is a snapshot of the various features loaded into every Safe 'n' Sure product:

BRAND	Pre-Loaded Cardboard Inserts	E-Z Glide Tab	Butterfly Seam	One Side Black - One Side See-Thru	Black On Both Sides
SAFE 'n' SURE OPT	YES	YES	YES	YES	NO
SAFE 'n' SURE DELUXE	NO	YES	YES	YES	NO
SAFE 'n' SURE ECONO	NO	NO	NO	NO	YES
SAFE 'n' SURE 3V	NO	NO	NO	YES	NO

Safe 'n' Sure DELUXE™ PSP Envelopes




Deluxe features our patented EZ-Glide tab, which makes loading and sealing quick and hassle free.

One side is semi-transparent so you can easily verify correct plate placement.

The easy-tear center seam lets you quickly and easily remove the plate, and you can load the plate into the reader without having to touch it too!

Deluxe comes 300 pcs per box.

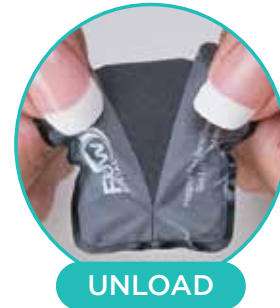

Safe 'n' Sure Deluxe

Size 0: 80101 Size 1: 80111 Size 2: 80121 Size 3: 80131
300/box

NEW! Safe 'n' Sure Deluxe in a Patented Wall Dispenser!

How do you make the best envelope even better? You offer it in a convenient wall dispenser!

New Safe'n'Sure Deluxe Wall Dispenser is a single use dispenser that can be mounted either on your plate processor or on a nearby wall. Easily dispenses one envelope at a time, so you can bag your plates just as they come out of your plate reader.

The dispenser comes with everything you need, and once finished, you simply order a new box.

Safe 'n' Sure Deluxe Wall Dispenser

Item: 80121-DBX Size 2 300/dispenser



Safe 'n' Sure Econo™

Econo features double sided black covering for complete protection against light exposure, plus our unique easy tear seam Envelopes load horizontally. After tearing open the envelope you have a thumb notch area to help remove the plate.

Safe 'n' Sure Econo

Size 0: 80100 Size 1: 80110 Size 2: 80120 Size 3: 80130

300/box



LOAD



SEAL



UNLOAD



Safe 'n' Sure Econo™ In a Wall Dispenser!

We've made our Econo even better by offering it in this handy roll -type dispenser. Simply mount the dispenser on a nearby wall or onto the plate reader using the adhesive tab provided. Then dispense as many envelopes as you need! Envelopes come on a roll for multiple dispensing.

Once done, simply re-order and replace your box /dispenser.

Safe 'n' Sure Econo Wall Dispenser

80120-DBX Size 2 300/box

Safe 'n' Sure 3V™

New Safe 'n' Sure 3V phosphor plate envelopes combine Safe 'N' Sure quality and a great low price!

These three "V"s stand for it's Value price, Vertical opening, and clear Viewing window. - making new Safe 'n' Sure 3V your best choice for an economically priced PSP envelope.

Safe 'n' Sure 3V

Size 0: 80122 Size 1: 80123 Size 2: 80124

300/box



LOAD




SEAL



UNLOAD



PSP Positioning Accessories

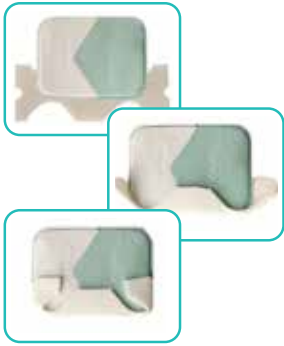

Deluxe Cushies™

These soft foam cushions are great for PSP envelopes, sensors, or film. Enhance patient comfort! Simply wrap the Deluxe Cushie around the image receptor. Its unique shape lets you cover all corners with a single Cushie. The V shaped notch helps to easily position the Cushie around corners. Sold in a dispenser of 300 pcs.

Deluxe Cushies™ 15810 300/box



SUPAs™

SUPA foam bite blocks hold your PSP envelopes at a 90° angle. For perfect alignment, use SUPA with our RAPID Positioning Rings and Arms.

SUPAs 16087 100/box



Sticky Bites™

Sticky Bites are our most comfortable bite block! Made of tasteless polyethylene foam, Sticky Bites easily affix to film or plate to help you capture the perfect image. Disposable and 100% Latex free.

Sticky Bites 15780 500/box



Stickems™

Adhesive bitewing tabs for use with either PSP barrier envelopes or film. For anterior, posterior and bite wings!

Stickems™ 15720 500/box



Cushies™

Cushies are foam tubes that create a soft edge for added patient comfort. Pre-slit to easily slip onto the long or short sides of the barrier envelope or film edge. Great for lower anteriors!

Cushies 15800 50/box



Loop EZ™

Made of heavy duty paper board, Loop-EZ are made to accommodate adult or child size films or PSP envelopes. Great for taking bite wings!

Loop-EZ	Adult Econo Pack	15725	1,500/box
	Adult:	15700	500/box
	Child:	15710	500/box

SMART™ Positioner

Flow's patented SMART system makes using phosphor plates just as easy, familiar and convenient as film.

Simply load your phosphor plate into the SMART Cassette and snap to close. The reusable cassette keeps phosphor plates safe from damage and protects them from moisture and light exposure. Place the cassette into the Universal Bite Block, attach an arm and an aiming ring, then position in the mouth.

SMART ensures perfect plate/cone alignment for every exposure. When you're done, take the plate out, process it and autoclave the cassette. The SMART system, including the cassettes, bite blocks, arms and rings, is completely autoclavable.

The SMART Kit 40854

Includes:

- 1 Anterior Aiming Arm
- 2 Anterior Aiming Rings
- 2 Posterior Aiming Arms
- 2 Posterior Aiming Rings
- 8 Universal Bite Blocks
- 20 Size 1 Cassettes
- 50 Size 2 Cassettes

Reorder Number

- 408570 (1 ea)
- 408590 (1 ea)
- 408580 (1 ea)
- 408600 (1 ea)
- 408700 (12 ea)
- 408500 (100 ea)
- 408520 (100 ea)

Additional parts can be ordered separately.



LOAD



SNAP CASSETTE IN BITE-BLOCK & EXPOSE



PROCESS



Clean Ups™

Perfect for cleaning your PSP plates, Clean Up's come individually wrapped. Single use wipes will remove fingerprints, dirt and other artifacts from your PSP plates.

Clean Ups 70107 50/box

RAPiD™ Positioner

The first color-coded PSP / film positioning system.

RAPiD bite blocks feature rounded corners and flex for easy loading and unloading. Aiming Rings and Bite blocks are color coded to make assembly faster and always accurate.

The RAPiD positioning system is available in a variety of kits, and all parts are sold separately.

For hand-held X-Rays like the Nomad* we make short versions of our RAPiD positioning arms. All kits are available with the shorty arms simply by adding a "-S" to our item numbers.



RAPiD Kit with Bite Wing 408560

RAPiD Shorty arms are great for hand-held x-rays!



RAPiD Kits

All Kits available with short positioning arms. Add a "-S" to the end of the item numbers below if you want to order a kit with short arms. Example: 408560-S

RAPiD Paralleling Kit with Bite-Wings 408560

- 4 Aiming Rings: 2 Anterior/Bite Wing, 2 Posterior
- 4 Indicator Arms: 1 Anterior, 2 Posterior, 1 Bite Wing
- 16 Bite Blocks: 4 Anterior, 8 Posterior, 4 Bite Wings

RAPiD Paralleling Kit without Bite Wings 409280

- 3 Aiming Rings: 1 Anterior, 2 Posterior
- 3 Indicator Arms: 1 Anterior, 2 Posterior
- 18 Bite Blocks: 6 Anterior, 12 Posterior

RAPiD Bite Wing Kit with #2 Horizontal Bite Blocks 409260

- 1 Aiming Ring, 1 Indicator Arm, 4 Bite Blocks

*Note: You cannot chemically or dry heat sterilize or use autoclave drying cycle.

Anterior Part Numbers

1 Ant. & Bite Wing Ring	408590
25 Anterior Bite Blocks	408610
1 Anterior Arm	408570
1 Shorty Arm	408570-S



Posterior Part Numbers

1 Posterior Aiming Ring	408600
25 Posterior Bite Blocks	408620
1 Posterior Arm	408580
1 Shorty Arm	408580-S



Anterior & Posterior

Bite Blocks 10 of each	408630
------------------------	--------

Bite Wing Part Numbers

6 #2 Hor.	409290
6 #2 Vert.	420410
6 #0/1 Hor.	409330
6 #3 Hor.	409310
1 Bite Wing Arm	409270
1 Shorty Arm	409270-S



Endodontics

Endodontic X-Ray Positioner Kit	40303
1 Aiming Ring	40310
1 Indicator Arm	40311
3 Bite Block	40312



Bite Block Baggies™

Why autoclave your Positioner parts or CBCT bite blocks when you can bag them instead? Bite Block Baggies are made to slide over the bite block and positioning bar. Prevent cross contamination the easy, inexpensive way. After each use simply re-bag. Economically priced, Bite Block Baggies come in a handy self-dispenser box. Each bag is 80 x 120 mm (3.14" x 4.72")



Bite Block Baggies 70104 500/box



Perfect Fit™ Adjustable Camera Sleeve

Quickly and easily achieve a tight, custom-fit on virtually any size camera head!

The only sleeve you'll ever need! Quick and easy! Your image is always crystal clear and sharp!

Perfect Fit features two-way customizable tabs. Cinch tight horizontally, cinch tight vertically, wrap the adhesive tab around and seal. With Perfect Fit there's no excess material to throw away. Waste is reduced and, at less than 17¢ a sleeve, your savings are increased!



Perfect Fit 80125 150 sleeves



Slide On



Fold Excess



Wrap & Stick

Comfee-C™ Composite Dispenser Sleeves

New Comfee-C sleeves are a perfect fit for your composite dispensers. Simply slide over the barrel and you're ready to use! Our unique design makes it easy to remove and replace after use.

Economically priced, our Comfee -C composite dispenser sleeves come in an easy to use dispenser box.



Comfee-C Box 80300 100 sleeves/box



Comfee Bytes™ Single Use Bite Covers

NEW! CBCT Users: Be kind to your patients and the environment with new COMFEE BYTES!


These Single use Bite Covers are a comfortable and inexpensive replacement for the Sirona* XIOS foam bite blocks, Sirona item 6141449. May also be used with Morita Verview, Planmeca Promax 3D, Sirona's XG3D, Axeos, and Vatech's Pax-i-3D.

Easy to use and completely biodegradable, new Comfee Bytes are the environmentally friendly and economical alternative to Sirona's* XiosFoam bite blocks. Comfee Bytes fold over the fixed bite block on your CBCT to form a soft, comfortable bite plate for patients to bite on.

Comfee Bytes will work with almost any 3D CT machine that features a bite block for patient positioning.



Comfee Bytes 81499 100/box

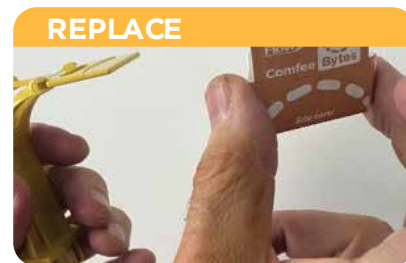

For Bite Block Baggies see page 7. For CBCT aprons see pages 14-15

Comfee Ears™ Disposable Ear Sleeves

Works with any CBCT X-Ray that includes ear positioners. Simply slide over your ear positioners and replace after every patient. Give your patient the peace of mind they deserve.

Comfee Ears 81500 500/box





Snuggies™ & Snuggies™ 2.0 Custom-Fit Sensor Sleeves

Protect your sensors and increase patient comfort with SNUGGIES adjustable sensor sleeves.

Snuggies work with any size sensors. Snuggies 2.0 is preferable for larger size sensors.

Simply slide your sensor inside the Snuggies sleeve, peel off the cover to the adhesive tab and wrap tight. No excess material for your patients to gag on. Your sleeve is always firmly in place over the sensor.

Regular Snuggies	80222	500/box
Snuggies 2.0 for Larger Sensors	80223	500/box

Deluxe Comfees™ Sensor Sleeves

Deluxe Comfees feature a super soft vinyl sleeve that makes loading sensor easy and effortless. No more struggling to get the sensor inside the sleeve. The added strength of the Deluxe Comfee sleeve protects you and your patients from accidental ripping.

Rounded corners for added patient comfort too. Tear-Away seams makes removing sensors quick and effortless.

Deluxe Comfees come in two sizes; Large (fits most size 2 sensors) and Medium (fits most size 0 and size 1 sensors). Each available in boxes of 100 and 500 pieces.

Comfees Deluxe	Large: 80213	Medium: 80212	100/box
Comfees Deluxe	Large: 80210	Medium: 80200	500/box



Textured Vinyl

Econo Comfees™ Sensor Sleeves

Econo Comfees are thinner, more form fitting sleeves than our Deluxe. Econo are economically priced too!

Made of LDPE material, Econo's feature tear-away seams to make removing your sensor quick and effortless.

Comfees Econo	Large: 80211	Medium: 80201	500/box
----------------------	--------------	---------------	---------



Smooth Seams

Digi Loops™

Digi Loops™ are made of heavy duty paperboard that simply 'loop' around the sensor. They're color-coded so you can easily pick out the proper size. Large and Medium hold most size 2 sensors. Small is suitable for most size 1 and size 0 sensors. Packaged 300 per box. Also available in an Assortment with 100 of each size.

Digi Loops	300/box		
Large: 15732	Medium: 15731	Small: 15730	Assortment 15733
Yellow color	Green color	Blue color	





Sensibles™ - Sensor Positioning System



Locking bumpers to hold bite block in multiple locations!

*Patent Pending

New and Improved design! Works with any size sensor! New softer corners for a comfortable patient experience!

Our patented Ratcheting bite block locks down to form a custom-fit around any size sensor in either the horizontal or vertical position. To remove the sensor, simply press the quick release button.

Sensibles also feature a floating bite plane so that you can move the sensor to the mid-line of your x-ray or make room for a shallow palette. The only sensor holder that can take Anteriors, Posteriors and Bite wings all at the same time!

Sensibles come in 2 sizes; large for most size 1 and 2 sensors in horizontal and vertical positions and Medium for most 0, 1, and 2 sensors in the horizontal position. Available in starter kits or each part may be purchased individually.

All parts autoclavable.



Using Sensibles is Easy As 1-2-3:



1 Attach Positioning Bar

2 Add Aiming Ring

3 Snap In Sensor

Sensibles Deluxe w/Bite Wing

- 6 Large & 10 Medium Bite Blocks
 - 2 Anterior Positioning Arms
 - 2 Anterior/Bite Wing Aiming Rings
 - 2 Bite Wing Positioning Arms
 - 2 Posterior Positioning Arms
 - 2 Posterior Aiming Rings
 - 50 Deluxe Comfees Samples
- 40915

Sensibles Deluxe Starter Kit

- 4 Large & 8 Medium Bite Blocks
 - 2 Anterior Positioning Arms
 - 2 Anterior Aiming Rings
 - 2 Posterior Positioning Arms
 - 2 Posterior Aiming Rings
 - 50 Deluxe Comfees Samples
- 40900

Sensibles Econo Starter Kit

- 2 Large & 2 Medium Bite Blocks
 - 1 Anterior Positioning Arm
 - 1 Anterior Aiming Ring
 - 1 Posterior Positioning Arm
 - 1 Posterior Aiming Ring
 - 25 Deluxe Comfees Samples
- 40901

Sensibles Refill Kit 40910

- 12 Large (Two-Tone) Bite Blocks

Sensibles Refill Kit 40909

- 12 Medium Bite Blocks

1 Anterior Positioning Arm 408570

1 Anterior/Bite Wing Aiming Ring 408590

1 Posterior Positioning Arm 408580

1 Posterior Aiming Ring 408600

1 Bite Wing Positioning Arm 409150



Cable Saver™

Cable Saver™ protects your sensor cables from being stepped on, rolled over, twisted, and tangled. Easy to install and easy to keep clean, Cable Savers will add years to the life of your sensor cables. Available in 8' lengths. Simply cut to size!

Cable Saver 80500 8' length



Uni-Verse-All™

Choose the adjustable bite block that best fits your sensor, snap into place on the Uni-Verse-All positioner, and you're good to go! Comes with aiming rings too. The aiming rings rotate into one of several positions to align your ring to the sensor in all Anterior, Posterior or bite wing positions!

Lightweight for added patient comfort. All parts autoclavable. Available in a starter kit. All parts sold separately too.

Uni-Verse-All Kit 40916

Kit contents:

- 2 Uni-Verse-All positioners
 - 4 Medium size Sensible's sensor holders
 - 4 Large size Sensible's sensor holder
 - 2 Aiming Rings
- (1) Instrument Only 40917
 - (12) Large Sensor holders 40910-U
 - (12) Med. Sensor Holders 40909-U
 - (1) Aiming Ring (purple) 408590



Bullseye™

Aiming ring for hand-held X-rays

Users of hand-held X-Rays like the Nomad* know how difficult it can be to maintain alignment between the positioner and the hand-held X-Ray. Bullseye is your answer! Simply peel and stick your Bullseye aiming ring onto the clear acrylic shield on your handheld X-Ray. Then the positioning arm of your sensor positioner passively nests into the Bullseye to complete your alignment.

Re-Usable. Economical. Works with all portable X-Rays that include an acrylic shield.


Bullseye Aiming Rings 408800 *Single Unit*



Peel off adhesive backing



Stick onto the acrylic ring

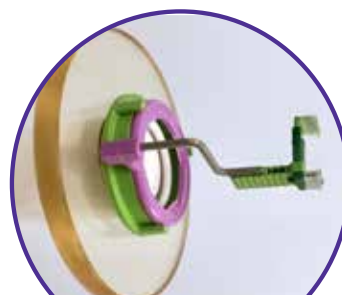

Nest your aiming ring and X-Ray!

All-Bite™

A one-size fits all sensor bite block for use without an aiming ring. Simple slide your sensor into the All Bite and snap the All Bite into the bite block. The bite block snaps in to form either a horizontal or vertical bite wing. Single use.

All-Bite 409155 30/box



Bite block attached for Horizontal Bite Wings



Bite block attached for Vertical Bite Wings





Dosimeter Badge Service

Keep you and your staff safe

TLD is state-of-the-art technology. Four elements are incorporated into our Dosimeters, so we measure down through three layers of tissue. And our processing software eliminates any guesswork or errors.

Here's how to keep you and your staff safe:

Everyone receives his or her own individual badge, to be worn at all times while at work. These lightweight badges measure 1" x 2" and clip easily onto clothing. Each badge is assigned a unique identification code.

Every 3 months a new replacement badge arrives in the mail, along with a return envelope. You send us back your old badge. We generate a detailed report and it's sent to you for your records. And if the readings ever show exposure in excess of published NRC safety and standards, your Flow service staff will contact you immediately by phone!

In addition to personal safety concerns, OSHA and many state agencies require that personnel use x-ray monitoring devices when operating, or in proximity to, any radiation-emitting equipment. Failure to comply may result in fines. And having individual documentation on the history of your staff's exposure to x-rays may prove vital in case of litigation any time in the future.

Enrollment is quick and easy:

- 1) Select who at your office will be the main contact and have them fill out the form below. If you need to list more names, copy the form.
- 2) Fax the form to (631) 242-1001
- 3) You will receive a return fax confirmation within 72 hours. Your first delivery of badges will arrive within 10 business days.

Your subscription includes automatic badge replacement every 3 months, for one year, reply envelopes for returning badges for processing, and copies of our quarterly reports (including reprints at no charge).



X-RAY BADGE SERVICE ORDER FORM

Refer to Stock Number 16102

Fax this form to 631-242-1001

Acct Number: _____ Clinic: _____ Dr.: _____

Address: _____ City: _____ State: _____ Zip: _____

Contact: _____ Phone: _____ Email: _____

Yes, I want to switch my service. My contract ends: ___/___/___ Yes, I want Automatic Annual Renewal

1. _____

Employee ID Number: _____

2. _____

Employee ID Number: _____

3. _____

Employee ID Number: _____

4. _____

Employee ID Number: _____

5. _____

Employee ID Number: _____

6. _____

Employee ID Number: _____

Exposure™ Articulating Film

Exposure's built-in handle eliminates the need for forceps. Easily exposes high occlusal contacts on both arches. Use for one tooth or an entire quadrant, on wet or dry teeth, even metal and non-metallic restorations!

Super thin Exposure Articulating Film makes it easy to identify and adjust high contacts. Each film measures 4 1/2 x 2 3/4".



Exposure Articulating Film 24300 200/box



Safe-T-First™ Face Shield

Fog-free, latex-free, one-size-fits-all full face shield

Safe-T-First face shields provide full face and eye protection. Unlike other face shields ours is fully adjustable for maximum safety and comfort. And they're lightweight so they can be comfortably worn for hours. The viewing area is both anti-fog and distortion free. Each shield come with a forehead cushion, (latex free), for maximum comfort. Economically priced, Safe-T-First face shields are individually wrapped.

Safe-T-First Face Shields WFS-01 24/box



Green, Lead-Free Protective Aprons

Flow's LEAD-FREE aprons make it easy to be green. Safe for the environment, lighter than conventional lead aprons! Made from materials that mimic the properties of lead, our lead-free garments are available in 6 aprons styles as well as a detachable thyroid collar, and come in 4 scintillating patterns.

There's also an antimicrobial cover that helps prevent cross contamination.

Flow lead free aprons are made to a 0.25 mm lead equivalency, and are 25% lighter than traditional leaded aprons. Aprons include a tear-resistant anti-slide backing and come with hanging loops to extend the life of your apron.



Separate Thyroid Collar in Paw Prints



CBCT Pano Cape without Collar


Panoramic Apron in Shimmering Green



How to order your Flow No-Lead Apron

Select your Apron style and cover pattern or color. Use the Stock number shown on the opposite page and add your cover fabric number. Then add a -NL to denote "no lead". **Example:** If you want an Adult Bib Apron in Shimmering Green, you would order item 75050NL-38.

Description	Size	Stock Number
Adult Bib Apron without Collar	24" x 26"	75050NL (+ color #)
Adult Bib Apron with Attached Collar	24" x 26"	75051NL (+ color #)
Child Size Bib Apron with Attached Collar	20" x 20"	75053NL (+ color #)
CBCT Pano Cape with Collar	23.5" x 23"	75090NL (+ color #)
CBCT Pano Cape without Collar	24" x 22.5"	75091NL (+ color #)
Adult Panoramic Apron	24" x 27"	75061NL (+ color #)
Child Panoramic Apron	20" x 20"	75062NL (+ color #)
Thyroid Collar	5" x 27"	75056NL (+ color #)
Apron Uni-Hanger		16401

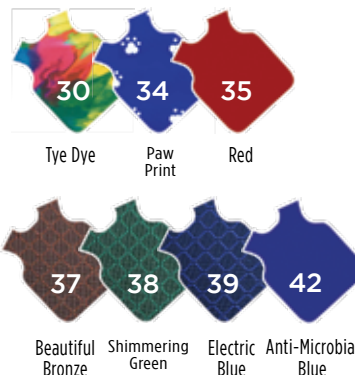



Lead Protective Aprons

Flexible lead protection never looked this good!

Flow's lead X-Ray aprons are the industry standard. Manufactured from flexible lead rubber and expertly sewn between an attractive easy to clean vinyl cover and a tear-resistant non-slide backing, our aprons combine great quality and economical prices.

Our lead aprons provide 0.3mm protection, come in 8 colors and 6 styles. Thyroid collars can be bought separately or sewn into the apron.

How to order your Flow Lead Apron

To order, select from any of the colors shown on the left. Use the stock number below and add a fabric number after it. **Example:** If you want a Bib Apron, without a thyroid collar, in the Teal fabric, the correct product number is 75050-14.



Description	Size	Stock Number
Adult Bib Apron without Collar	24" x 26"	75050 (+ color #)
Adult Bib Apron with Attached Collar	24" x 26"	75051 (+ color #)
Child Size Bib Apron with Attached Collar	20" x 20"	75053 (+ color #)
CBCT Pano Cape with Collar	23.5" x 23"	75090 (+ color #)
CBCT Pano Cape without Collar	24" x 22.5"	75091 (+ color #)
Adult Panoramic Apron	24" x 27"	75061 (+ color #)
Child Panoramic Apron	20" x 20"	75062 (+ color #)
Thyroid Collar	5" x 27"	75056 (+ color #)
Apron Uni-Hanger		16401

Tranzform Adjustable Apron Hanger



Improperly store your X-Ray aprons and you can get cracks in the protective material that changes "protective" into a guessing game. New Tranzform Adjustable Apron Hanger changes all that. Tranzform is the first apron hanger that adjusts to the size and shape of your apron.

And you can easily store 2 Aprons on a single Tranzform Hanger.

Tranzform Apron Hanger 16220 1 per box

Uni-Hanger



Strong all metal construction. Holds a single apron

Uni-Hanger Apron Hanger 16301 1 per box



Silver D: D Speed X-Ray Film

Flow Silver D is a high contrast film especially designed for today's most discerning dentist. This higher contrast delivers greater resolution and detail.

Now That's an Image!

In independent test results, Silver D delivered as much as 21% greater contrast and over 40% better resolution than other competitor's films.

What's more, our unique formulation stands up to oxidized chemistry better than other films. That means you'll get more films out of each gallon of developer used.

That's what on the inside. Outside, our vinyl packets are thinner, softer, and more flexible. Your patients will love our softer edges, your assistants will enjoy the added flexibility and easy-opening tabs!

Silver D has the same exposure settings as other D speed films and can be processed manually or automatically.



Econo Pak
600/Box

Value Pak
150/Box

Fresh Pak
100/Box

*Ultra-Speed and Insight are registered trademarks of the Carestream Dental/Eastman Kodak Company.

Size 2

1 1/4" x 1 5/8"
31 x 41 mm

DV-57 DV-58
FV-57 FV-58

Size 0

7/8" x 1 3/8"
22 x 35 mm

DV-53 DV-54
FV-53 FV-54

Speeds

Flow Dental's dental films are available in D speed (similar to Carestream Ultra-Speed*) or F speed (similar to Carestream Insight*) which requires roughly half the exposure time of D speed.

Specifications

Each Flow Intraoral film packet contains films sandwiched between black paper and lead foil. The film is embossed with a dot and the packets are color coded to indicate speed type and to ensure the correct side of the film is placed to the tube. Flow film is available in five intraoral and two extraoral sizes, and can be processed manually or automatically.



Flow Dental Silver D Film (D Speed) Exposure Time In Seconds

		Exposure Factors	60kV 7mA	65kV 7.5mA	65kV 8mA	70kV 7mA	70kV 8mA	80kV 10mA
REGION	Maxillary	Anterior	0.55	0.32	0.30	0.27	0.24	0.10
		Premolar	0.73	0.43	0.40	0.37	0.32	0.13
		Molar	0.82	0.48	0.45	0.41	0.36	0.14
		Bite Wing	0.55	0.32	0.30	0.27	0.24	0.10
	Mandibular	Anterior	0.46	0.27	0.25	0.23	0.20	0.08
		Premolar	0.50	0.29	0.28	0.25	0.22	0.09
		Molar	0.55	0.32	0.30	0.27	0.24	0.10
		Bite Wing	0.55	0.32	0.30	0.27	0.24	0.10

Exposure times shown are for an 8 inch Source-image distance. For a 12 inch cone, multiply the exposure time by 2.25. For a 16 inch cone, multiply the exposure time by 4. To convert times to impulses, multiply the exposure times by 60. NOTE: For large patients, increase time by approximately 25%. For children, reduce exposure time by approximately 30%. For edentulous patients, reduce exposure time by approximately 25%. For DC x-ray generators, reduce exposure time by about one-third.



Silver D: D Speed X-Ray Film



DV-58 Film

D Speed Size 2 Adult Film

DV-58 D-speed, size 2 intraoral film is packaged in soft, supple vinyl packets for greater patient comfort. Perfect for a broad range of x-ray exams, including anterior, posterior, and bite-wing exams. DV-58 is available in three package sizes.

Econo Pak®
18204

600 films per box

Value Pak®
18200

150 films per box

Fresh Pak®
18600

100 films per box



DV-57 Double Film

D Speed Size 2 Adult Film - 2 Film Packets

All the great qualities of Silver D Film in a double film packet. One exposure = Two films! Perfect for referral situations or when a film duplicate is required for insurance processing. Twin D has the same exposure settings as regular Silver D, and can be processed manually or automatically.

Value Pak®
18201

130 films per box

Fresh Pak®
18601

100 films per box



DV-54 & DV-53 Double Film

D Speed Size 0 Child Film

All the great image qualities of Flow's Silver D in a size 0, child film size. Pedo Pak is perfect for kids ages 3 through 8. Pedo Pak is available in single film packets and double film packets, (Double D Pedo Pak). Pedo Paks can be processed automatically or manually.

Pedo Econo Pak®
18206

400 single films per box

Pedo Pak®
18235

100 single films per box

Double D Pedo Pak®
18236

100 double films per box



Xpress: F Speed X-Ray Film

As fast as any intraoral film on the market, Flow F speed films combine reduced radiation exposure with superb contrast and image resolution.

Don't Sacrifice Image Quality for Value, Speed & Comfort!

Available in sizes 2 and 0, in single and double film packets, Flow F speed takes conventional dental radiography to a new level in terms of dose minimization.

Inside: A fast film that doesn't sacrifice image quality for reduced radiation exposure! Outside: softer, thinner, more flexible vinyl packets your patients will love! Flow F speed films require approximately half the exposure time of D speed (see the exposure table below) and can be processed manually or automatically. Flow F speed has the same exposure settings as other F speed films making it easy to switch.



*Ultra-Speed and Insight are registered trademarks of the Carestream Dental/Eastman Kodak Company.

Size 2
 1 1/4" x 1 5/8"
 31 x 41 mm
 DV-57 DV-58
 FV-57 FV-58

Size 0
 7/8" x 1 3/8"
 22 x 35 mm
 DV-53 DV-54
 FV-53 FV-54

Speeds

Flow Dental's dental films are available in D speed (similar to Carestream Ultra-Speed*) or F speed (similar to Carestream Insight*) which requires roughly half the exposure time of D speed.

Specifications

Each Flow Intraoral film packet contains films sandwiched between black paper and lead foil. The film is embossed with a dot and the packets are color coded to indicate speed type and to ensure the correct side of the film is placed to the tube. Flow film is available in five intraoral and two extraoral sizes, and can be processed manually or automatically.



Flow Dental Xpress Film (F Speed) Exposure Time In Seconds

		Exposure Factors	60kV 7mA	65kV 7.5mA	65kV 8mA	70kV 7mA	70kV 8mA	80kV 10mA
REGION	Maxillary	Anterior	0.25	0.14	0.14	0.12	0.11	0.04
		Premolar	0.33	0.19	0.18	0.16	0.14	0.06
		Molar	0.37	0.22	0.20	0.19	0.16	0.06
		Bite Wing	0.25	0.14	0.14	0.12	0.11	0.04
	Mandibular	Anterior	0.21	0.12	0.11	0.10	0.09	0.04
		Premolar	0.23	0.13	0.12	0.11	0.10	0.04
		Molar	0.25	0.14	0.14	0.12	0.11	0.04
		Bite Wing	0.25	0.14	0.14	0.12	0.11	0.04

Exposure times shown are for an 8 inch Source-image distance. For a 12 inch cone, multiply the exposure time by 2.25. For a 16 inch cone, multiply the exposure time by 4. To convert times to impulses, multiply the exposure times by 60. NOTE: For large patients, increase time by approximately 25%. For children, reduce exposure time by approximately 30%. For edentulous patients, reduce exposure time by approximately 25%. For DC x-ray generators, reduce exposure time by about one-third.



Xpress: F Speed X-Ray Film



FV-58 Film

F Speed Size 2 Adult Film

Flow's FV-58 is an F speed, size 2, intraoral dental film in soft vinyl packets. FV-58 requires approximately half the exposure time of DV-58.

FV-58 can be used for anterior, posterior and bite-wing x-ray exams in adults, and posterior and bite-wing exams in children 9 years or older and can be processed manually or automatically.

Econo Pak®
18215
600 films per box

Value Pak®
18310
150 films per box

Fresh Pak®
18620
100 films per box



FV-57 Double Film

F Speed Size 2 Adult Film - 2 Film Packets

All the great image qualities of our super fast F speed Adult dental x-ray film in a double film packet. One exposure = two films!

Perfect for referral situations or when a duplicate film is required for insurance processing. Flow F Speed film uses half the exposure settings of our Silver D, and can be processed manually or automatically.

Value Pak®
18410
130 films per box

Fresh Pak®
18621
100 films per box



FV-54 & FV-53 Double Film

F Speed Size 0 Child Film

All the great image qualities of Flow's Xpress F in a size 0, child film size. Pedo Pak is perfect for kids ages 3 through 8. Pedo Pak is available in single film packets and double film packets, (Double F Pedo Pak). Pedo Paks can be processed automatically or manually.

Pedo Econo Pak®
18207
400 single films per box

Pedo Pak®
18635
100 single films per box

Double F Pedo Pak®
18636
100 double films per box



Butterfly Barrier Packet Films

Butterfly Barrier Packet Films are your preferred alternative to Carestream's Clinasept[®] film packets! Our easy-opening Butterfly packets come pre-loaded with your favorite Flow film, and protect you and the film from oral fluids. Enjoy the added safety and convenience while reducing your cost!

Butterfly's unique "wings" at the top of the envelope allow you to open the packet with minimum contact, keeping fingers free of blood and saliva. As an extra assurance of minimum contact, our envelopes open by pulling the "wings" apart from left-to-right as opposed to front-to-back. Once cleanly removed from its envelope, the remaining piece of film is a standard, easy-to-open, Flow film packet.

Butterfly films are available in the most popular sizes, styles and speeds: D and F Speeds, in Adult and Child sizes, and in single and double film packets.



	Speed	Films per Packet	Butterfly Name	Pieces per Box	Stock Number	Type of Film
Size 2	D-Speed	Single	DSS-58	100	19600	DV-58
		Double	DSS-57	100	19601	DV57
	F-Speed	Single	FSS-58	100	19620	FV-58
		Double	FSS-57	100	19610	FV-57
Size 0	D-Speed	Single	DSS-54	75	19235	DV-54
		Double	DSS-53	75	19236	DV-53
	F-Speed	Single	FSS-54	75	19635	FV-54
		Double	FSS-53	75	19636	FV-53

**Clinasept is a registered trademark of Carestream Dental/Eastman Kodak Company*

Bite-Chex Paper Bite Wing Film Size 3

Size 3 Bite wing film made in paper packets.

Bite-Chex Paper Bite Wing 18060
100/box



Single Chip Hangers

NEED DESCRIPTION

Single Chip Hangers 16210
12/box



Value Chem[®] Concentrated Developer and Fixer

Flow Developer has a Phenidone base for a higher concentration of developer, yielding better radiographs and slower oxidation of the solution. Flow Fixer is designed to fix a film in three minutes. A bottle of concentrate makes one gallon of working solution.

Developer & Fixer 24220 1 each





Panoramic Film & Accessories



Green Team® 512 & 612

Full speed, high contrast, orthochromatic (green sensitive) dental x-ray film available in 5 x 12" and 6 x 12" sheets.

Flow's Green Team® is ideally suited for all panoramic and TMJ imaging. The spectral sensitivity of Green Team is perfectly matched for use with all brands of green emitting, rare earth, intensifying screens. Green Team has wide surface area silver halide crystals to yield enhanced image sharpness.

Green Team® 512 & 612

Green Team 512 - 5 x 12" Green Sensitive Pano Film **18095G-F**

Green Team 612 - 6 x 12" Green Sensitive Pano Film **18096G-F**

50 films per box



True Blue® 512 & 612

Full speed, blue sensitive dental x-ray film available in 5 x 12" and 6 x 12" sheets.

Flow's True Blue® Panoramic film yields excellent detail with high contrast. The blue spectral sensitivity of True Blue is perfectly matched for use with all brands of calcium tungstate (blue emitting) intensifying screens.

True Blue® 512 & 612

True Blue 512 - 5 x 12" Blue Sensitive Pano Film **18095-F**

True Blue 612 - 6 x 12" Blue Sensitive Pano Film **18096-F**

50 films per box

Kopy Kat Duplicator

This duplicator, when used with duplicating film, will copy films up to a 6 x 12" size. Simply place the film to be copied onto the glass screen of the duplicator. Then place a sheet of duplicating film on top, close the cover and run. Then process your duplicated film in your film processor.



Kopy Kat Duplicator **21520**

Flexible Panoramic Film Cassette

Flexible and easy to load, our panoramic cassettes are made from heavy duty light proof vinyl. Can be used with any brand film.



Flexible Panoramic Film Cassette

Intended only for Panoramic Corp panoramic x-ray machines. **11532 1 per box**

Slimline Illuminator

Perfect for your desk or wall-hung, Slimline Illuminators are only 4 3/4" deep. Great for viewing intra oral or panoramic films, the Slimline has a 6" x 13 1/2" Viewing area.

Comes with a 6" line cord and keyholes in the back for wall mounting.



Slimline Illuminators **21502**

Angular Safelight

The Angular Safelight is ideal where a direct light is wanted. Easy to install via 2 keyholes, the Angular Safelight includes a GBX type filter, light bulb, and a 6 foot cord with a three prong plug.

Specifications: 10" H x 7 1/2" W x 5 1/2" D.



Angular 5 x 7" Safelight **19103**

Angular 5 x 7" with auto room light cut off **19203**

Circular 3" dia. **19104**



Film Positioning Accessories



Deluxe Cushies™

These soft foam cushions are great for PSP envelopes, sensors, or film. Enhance patient comfort! Simply wrap the Deluxe Cushie around the image receptor. Its unique shape lets you cover all corners with a single Cushie. The V shaped notch helps to easily position the Cushie around corners. Sold in a dispenser of 300 pcs.

Deluxe Cushies™ 15810 300/box



SUPAs™

SUPA foam bite blocks hold your PSP envelopes at a 90° angle. For perfect alignment, use SUPA with our RAPID Positioning Rings and Arms.

SUPAs 16087 100/box



Sticky Bites™

Sticky Bites are our most comfortable bite block! Made of tasteless polyethylene foam, Sticky Bites easily affix to film or plate to help you capture the perfect image. Disposable and 100% Latex free.

Sticky Bites 15780 500/box



Stickems™

Adhesive biting tabs for use with either PSP barrier envelopes or film. For anterior, posterior and bite wings!

Stickems™ 15720 500/box



Cushies™

Cushies are foam tubes that create a soft edge for added patient comfort. Pre-slit to easily slip onto the long or short sides of the barrier envelope or film edge. Great for lower anteriors!

Cushies 15800 50/box



Loop EZ™

Made of heavy duty paper board, Loop-EZ are made to accommodate adult or child size films or PSP envelopes. Great for taking bite wings!

Loop-EZ	Adult Econo Pack	15725	1,500/box
	Adult:	15700	500/box
	Child:	15710	500/box



RAPiD™ Positioner

The first color-coded PSP / film positioning system.

RAPiD bite blocks feature rounded corners and flex for easy loading and unloading. Aiming Rings and Bite blocks are color coded to make assembly faster and always accurate.


The RAPiD positioning system is available in a variety of kits, and all parts are sold separately.

For hand-held X-Rays like the Nomad* we make short versions of our RAPiD positioning arms. All kits are available with the shorty arms simply by adding a "-S" to our item numbers.



RAPiD Kit with Bite Wing 408560

RAPiD Shorty arms are great for hand-held x-rays!

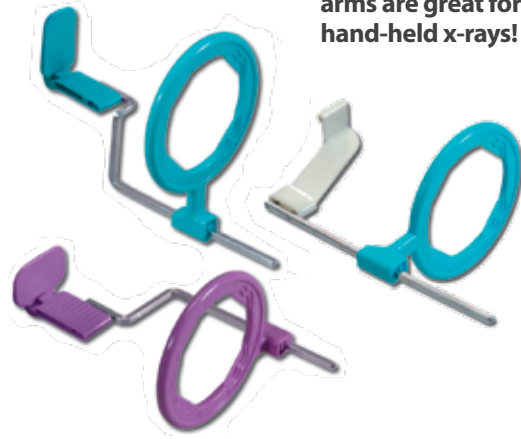

RAPiD Kits

All Kits available with short positioning arms. Add a "-S" to the end of the item numbers below if you want to order a kit with short arms. Example: 408560-S

RAPiD Paralleling Kit with Bite-Wings 408560

- 4 Aiming Rings: 2 Anterior/Bite Wing, 2 Posterior
- 4 Indicator Arms: 1 Anterior, 2 Posterior, 1 Bite Wing
- 16 Bite Blocks: 4 Anterior, 8 Posterior, 4 Bite Wings

RAPiD Paralleling Kit without Bite Wings 409280

- 3 Aiming Rings: 1 Anterior, 2 Posterior
- 3 Indicator Arms: 1 Anterior, 2 Posterior
- 18 Bite Blocks: 6 Anterior, 12 Posterior

RAPiD Bite Wing Kit with #2 Horizontal Bite Blocks 409260

- 1 Aiming Ring, 1 Indicator Arm, 4 Bite Blocks

*Note: You cannot chemiclave or dry heat sterilize or use autoclave drying cycle.



Anterior Part Numbers

1 Ant. & Bite Wing Ring	408590
25 Anterior Bite Blocks	408610
1 Anterior Arm	408570
1 Shorty Arm	408570-S



Posterior Part Numbers

1 Posterior Aiming Ring	408600
25 Posterior Bite Blocks	408620
1 Posterior Arm	408580
1 Shorty Arm	408580-S

Anterior & Posterior

Bite Blocks 10 of each	408630
------------------------	--------



Bite Wing Part Numbers

6 #2 Hor.	409290
6 #2 Vert.	420410
6 #0/1 Hor.	409330
6 #3 Hor.	409310
1 Bite Wing Arm	409270
1 Shorty Arm	409270-S



Endodontics

Endodontic X-Ray Positioner Kit	40303
1 Aiming Ring	40310
1 Indicator Arm	40311
3 Bite Block	40312

Bite Block Baggies™



Why autoclave your Positioner parts or CBCT bite blocks when you can bag them instead? Bite Block Baggies are made to slide over the bite block and positioning bar. Prevent cross contamination the easy, inexpensive way. After each use simply re-bag. Economically priced, Bite Block Baggies come in a handy self-dispenser box. Each bag is 80 x 120 mm (3.14" x 4.72")

Bite Block Baggies 70104 500/box



X-Ray Mounts

Flow Dental X-Ray mounts are manufactured of the finest archival materials. Flow makes four types of mounts, each available in popular styles. Whichever your preference, Flow mounts will give you decades of use.

Perfect Pocket™ Mounts

A great alternative to Rinn's E-Z View* masked mounts, Perfect Pocket mounts feature an innovative design loaded with benefits others don't have. Our material is thinner, making it easier to load films, and the vinyl is clearer than others, so you have a better view of your films. Every Perfect Pocket mount features rounded corners- no more finger cuts! With additional space for writing patient information.

Also available in Tear-Away - allowing you to use as many windows as you need.



Perfect Pocket Mounts

See chart on page 25 for item #s 100/box

Tear-Away Mounts

81-0500 80/box (960 windows)

Klear-Vue™ Cardboard Mounts

Three layers of cardboard form a recessed pocket and lip that hold film in place. Flow mounts are printed with room for patient data (name, date, notations). Mounts can be imprinted with personalized doctor information. Klear-Vue mounts are comparable to Rinn's EeZee mounts*.



Klear-Vue Cardboard Mounts

See chart on page 25 for item #s 100/box

Greylight™ Mounts

Similar to the Rinn EZ-Tab* and ADA* style mounts, our Greylight Mounts feature open window viewing. Films are held in place via embossed frames with raised borders that load easily and lock films firmly into place. Greylight mounts can be imprinted with custom doctor information.



Greylight Mounts

See chart on page 25 for item #s 100/box

Grip Right™ Plastic Mounts

Similar to the Rinn EeZee Mount* plastic and Greene* TransView Mounts. Grip Right Plastic mounts are made from a thin semi-rigid vinyl with die-cut tabs. Thinner than other mount styles, Grip Right mounts allow for compact storage of your patient records. Can be imprinted with custom doctor information.



Grip Right Plastic Mounts

See chart on page 25 for item #s 100/box

*ADA Mounts are a Registered Trademark of ADA Dental Products. *Rinn EeZee Mount, EZ-View, EZ Tab and Greene Transview are Registered Trademarks of Dentsply International. * Trollmount is a Registered Trademark of Trollplast.



X-Ray Mounts Reference Chart

Layout	Description		Perfect Pocket Mounts	Klear-Vue Cardboard Mounts	Grip Right Plastic Mounts	Graylight Plastic Mounts
	1H	1 Wind. Horiz. Size 2	81-6011	2500		
	2H	2 wind. Horiz. Size 2	81-6021	25203		602-203G
	2H-0	2 wind. Horiz. Size 0	81-4021	25608	23608	
	2H-BT	2 wind. Horiz. Size 3		25609	23206	
	3H	3 Wind. Horiz. Size 2	81-7021			
	4H	4 Wind. Horiz. Size 2	81-6040	25209	23209	604-209G
	4H-0	4 Wind. Horiz. Size 0		25609	23609	
	4H-2R	4 Wind. Horiz., 2 Rows, Size 2	81-6045		23239	6045-239G
	6H	6 Wind. Horiz., 2 Rows, Size 2	81-6061	25215		606-215G
	6-2V	6 Wind. 4 Horiz., 2 Vert. Size 2			23216	
	8H	8 Wind. Horiz., 2 Rows Size 2	81-6081	25219	23219	608-219G
	10B-2	10 wind B wing Size 2		25261		
	10BT-3	10 wind B wing Size 3				610-262G
	10-6V	10 Wind. 4 Horiz. 6 Vert. Size 2		25240		
	14-6V	14 Wind. 8 Horiz. 6 Vert. Size 2	81-8142	25225		614-225G
	14-6V-0	14 Wind. 8 Horiz. 6 Vert. Size 0				614B-625G
	16-6V	16 Wind. 10 Horiz. 6 Vert., Size 2		25226	23226	
	16-6V 2BT	16 Wind. 8 Horiz. 6 Vert. Size 2, B Wing		25227		
	18-6V	18 Wind. 12 Horiz. 6 Vert. Size 2	81-6180	25228	23228	618-228G
	20H	20 Wind. Horiz. Size 2	81-4567			618-230G
	20-12H "2" 8V-1	20 Wind. 12 Horiz., 8 Vert. Size 1				610-270G
	Roll Type	Tear-away Roll of 960 Pockets	81-0500			



133 Wolf Road, Battleboro, NC 27809

Phone: 800.356.9729

www.flowdental.com | info@flowdental.com